







Site-specific Physical Failure Analysis of 3D Systems using “Plasma” FIB

Remco Geurts, P.D. Carleson, R.J. Young, R. Routh, C. Rue, G. Franz, and L.F.Tz. Kwakman






FEI Company

EFUG 2011

Monday 3rd October, 2011

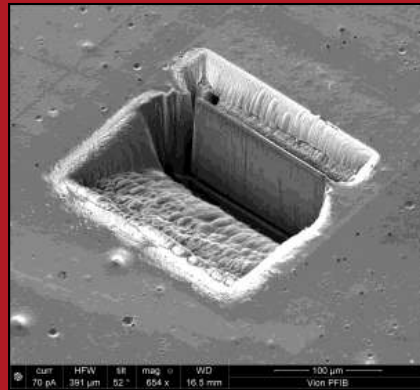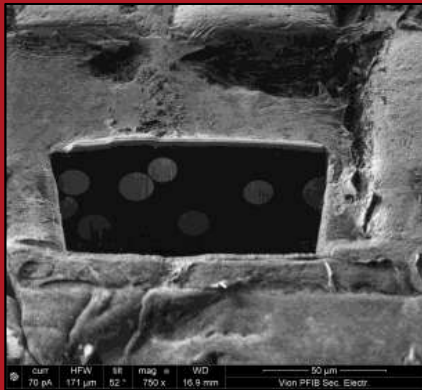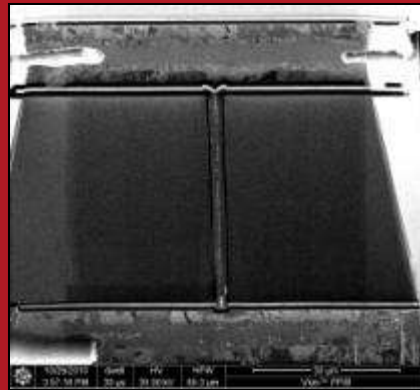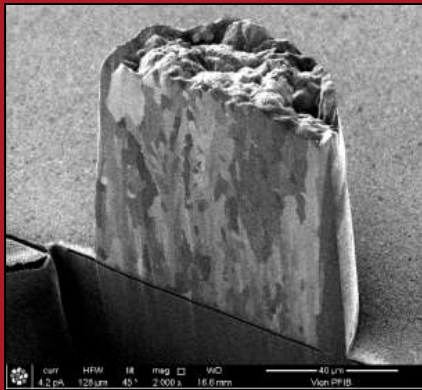


Plasma FIB Technology and Applications

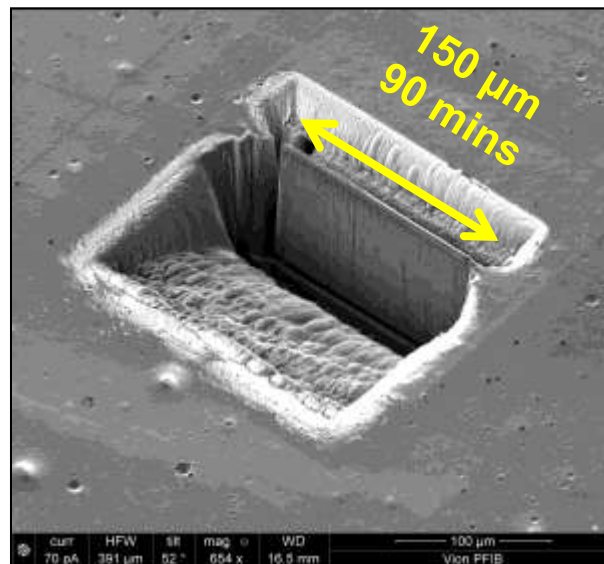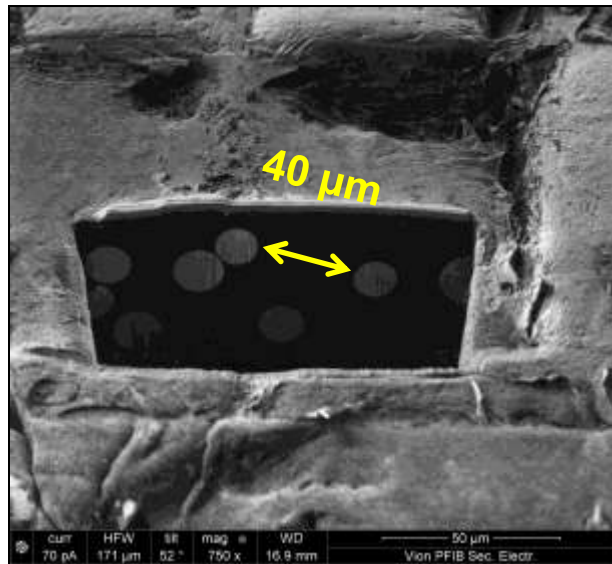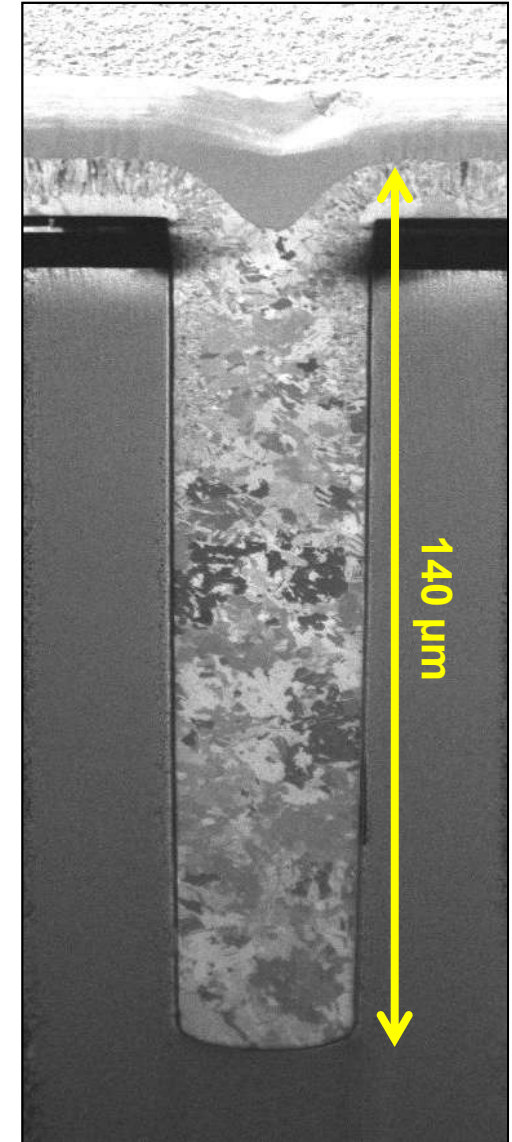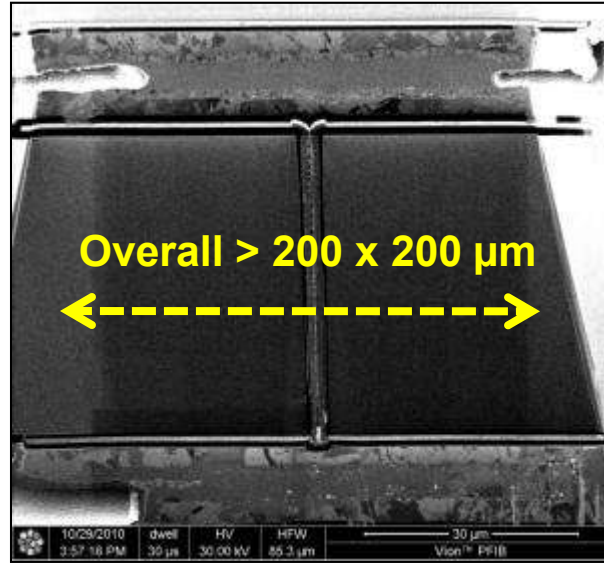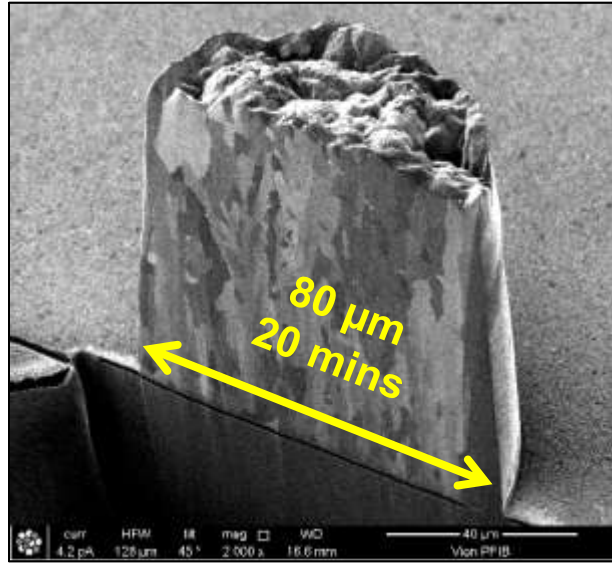

Focused ion beam sample preparation

Excellent for site-specific sample preparation


- Localized - positioning to nm level
 - Leave rest of device intact
 - Multiple locations on single device
- Any orientation
- No mechanical shock or tearing/smearing

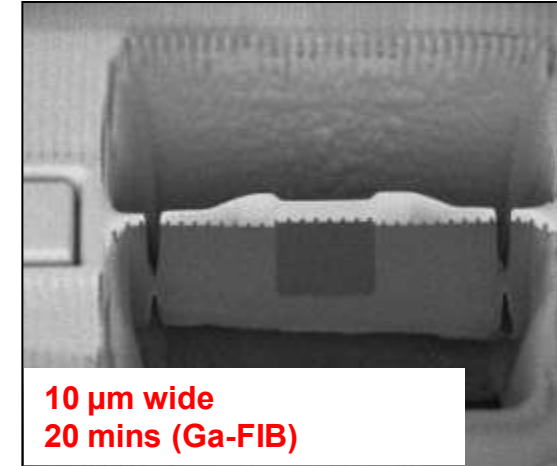

Generally uses gallium liquid metal ion source (LMIS)


Typical beam current range 1 pA to 20-65 nA

- $\sim 10^3 \mu\text{m}^3/\text{min}$ for silicon at 60 nA
- Short prep times for sections a few 10s of μm on a side
- 3D IC technology/packaging sections often $> 100 \mu\text{m}$
- Therefore, require new techniques focused on throughput and efficiency

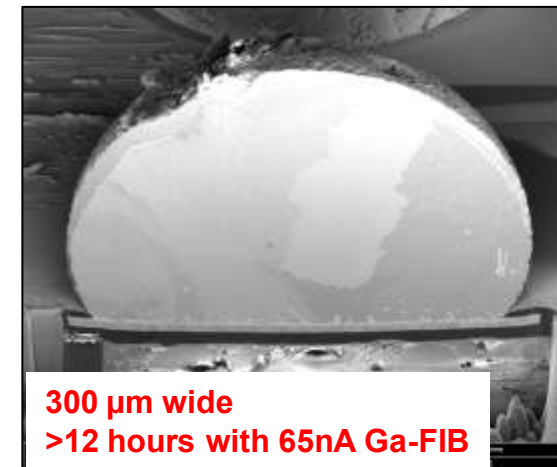


Why *Plasma* FIB: 20x faster than current FIBs

A system that provides unique and fast ion milling capabilities for rapid cross sectioning of features from 50 to 1000 microns.

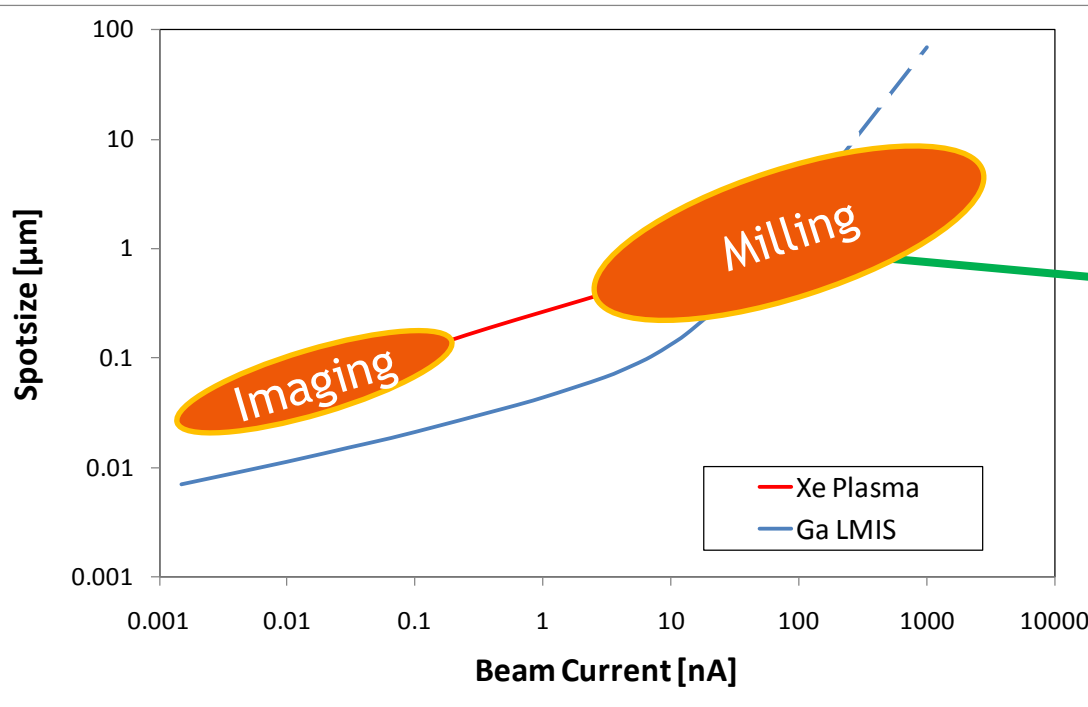

- ✓ High volume milling / high beam current
- ✓ Ga-FIB loses size advantage to plasma source as beam current goes above 50-60 nA
- ✓ Xe has high sputter yield, high brightness, and low energy spread
- ✓ No Ga contamination


Inductively coupled plasma (ICP) ion source

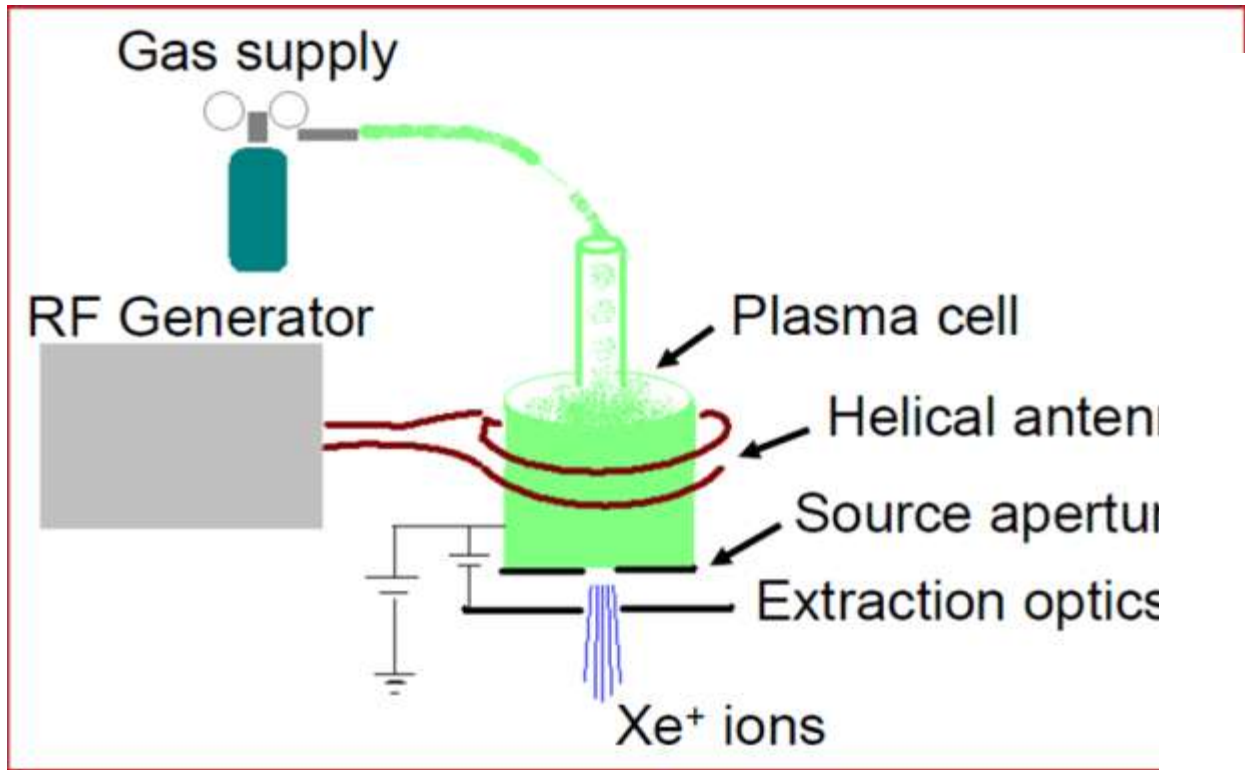

- Gas flows into plasma cell
- Helical antenna couples energy into plasma cell
- Electrons removed from atoms to form Xe⁺ ions
- Extraction optics accelerate ions into FIB column


Comparison of FIB sources

	Liquid Metal Ion Source (LMIS)	Inductively Coupled Plasma (ICP) Source
Descriptive Phrase	“Point source with low angular intensity”	“Broad source with high angular intensity”
Angular Intensity [mA sr⁻¹]	0.02	50
Virtual Source Diameter [nm]	50	15,000
Normalized Brightness [A sr⁻¹ m⁻² V⁻¹]	10 ⁶	10 ⁴
Optical Challenge	Above about 10 nA, spherical aberrations degrade FIB performance	Extreme demagnification required for small probe sizes (< 50 nm)

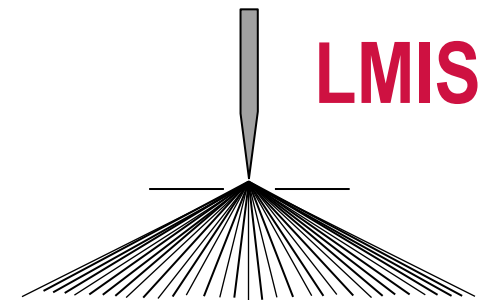

LMIS



Conical emission pattern

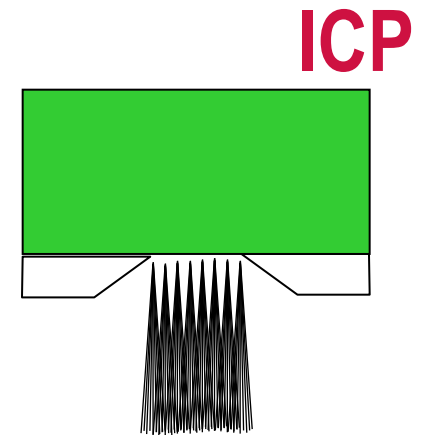

ICP

Collimated emission pattern


TSV- full section (depo, bulk, polish) in 20 mins

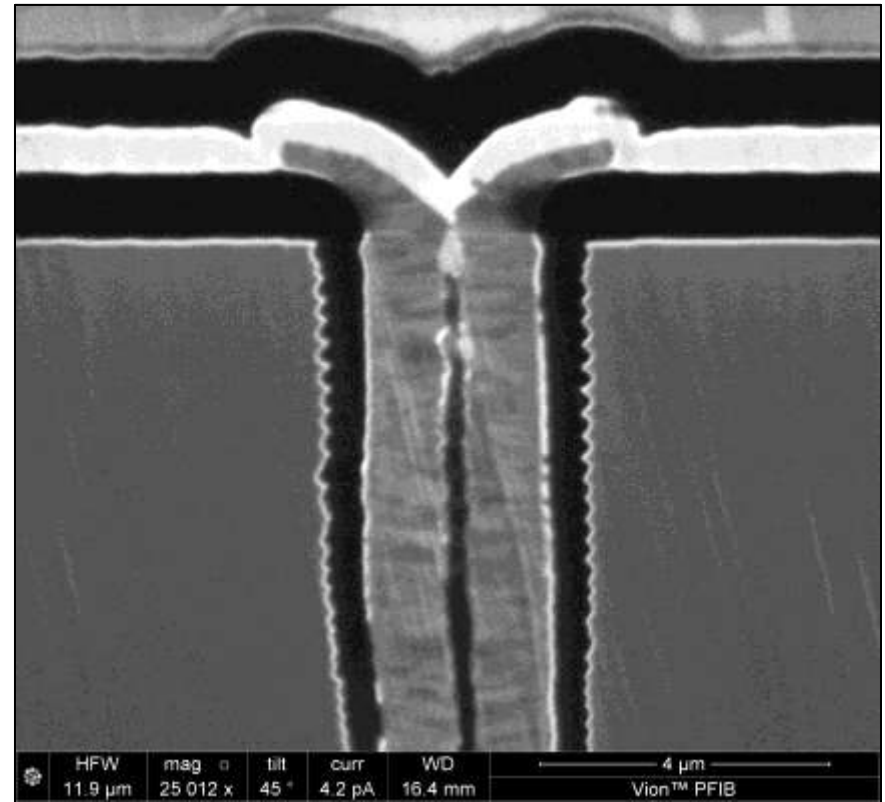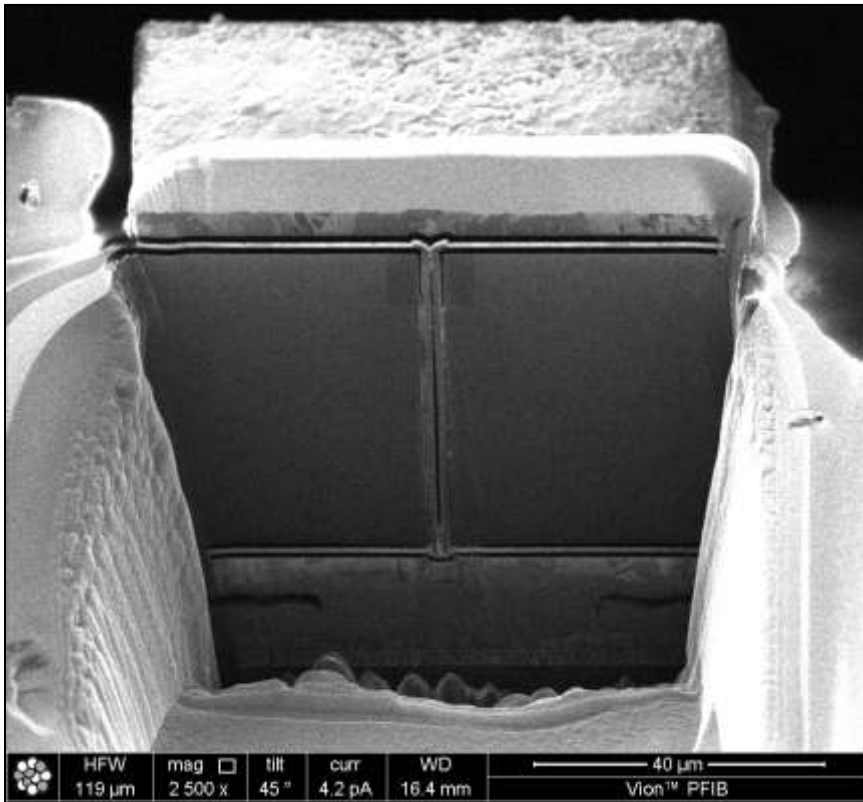

PFIB imaging provides the similar contrast mechanisms as standard Ga-FIB:



- Topographic contrast
- Materials contrast
- Passive voltage contrast
- Ion channeling contrast
- Smallest useful HFW: 10-15 μm
- Limited sputtering/image at low current (compared to using HFW of 2-3 μm with Ga)

Cross sectioning three-die stack

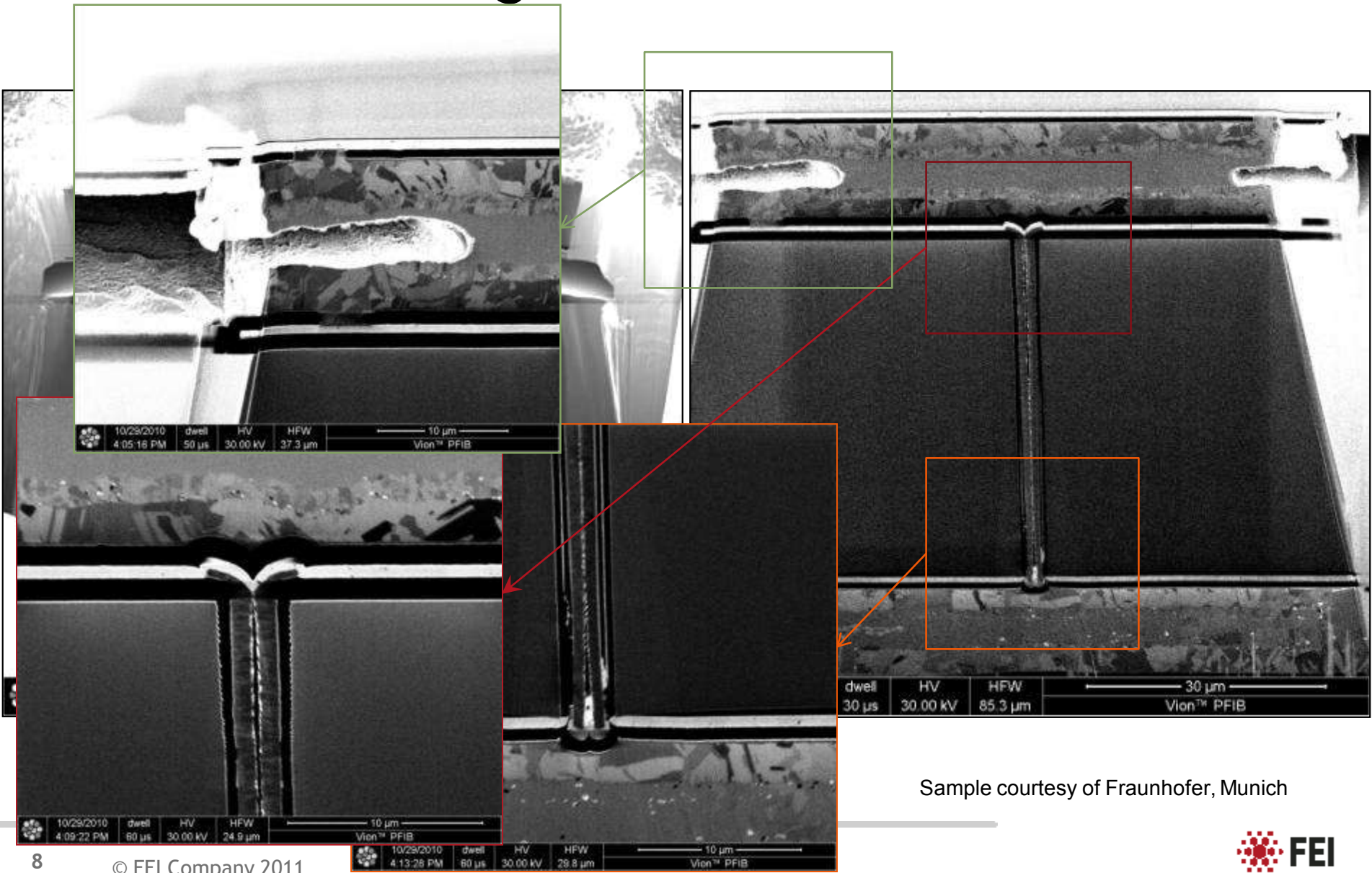


Sample courtesy of Fraunhofer, Munich

Details of SLID metallurgy (Cu-Sn-Cu)

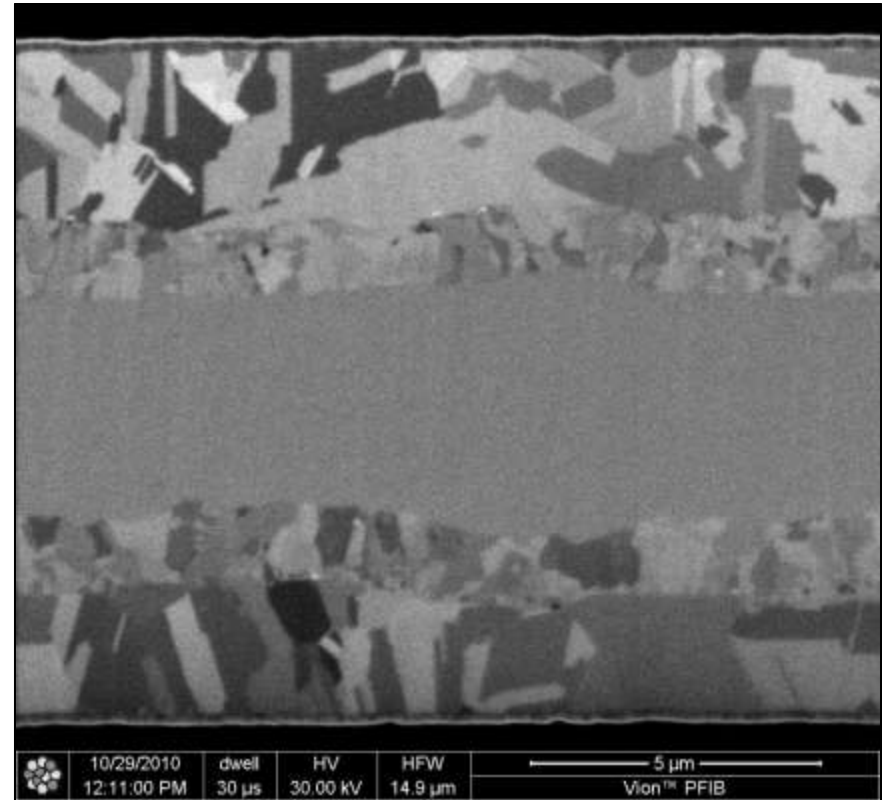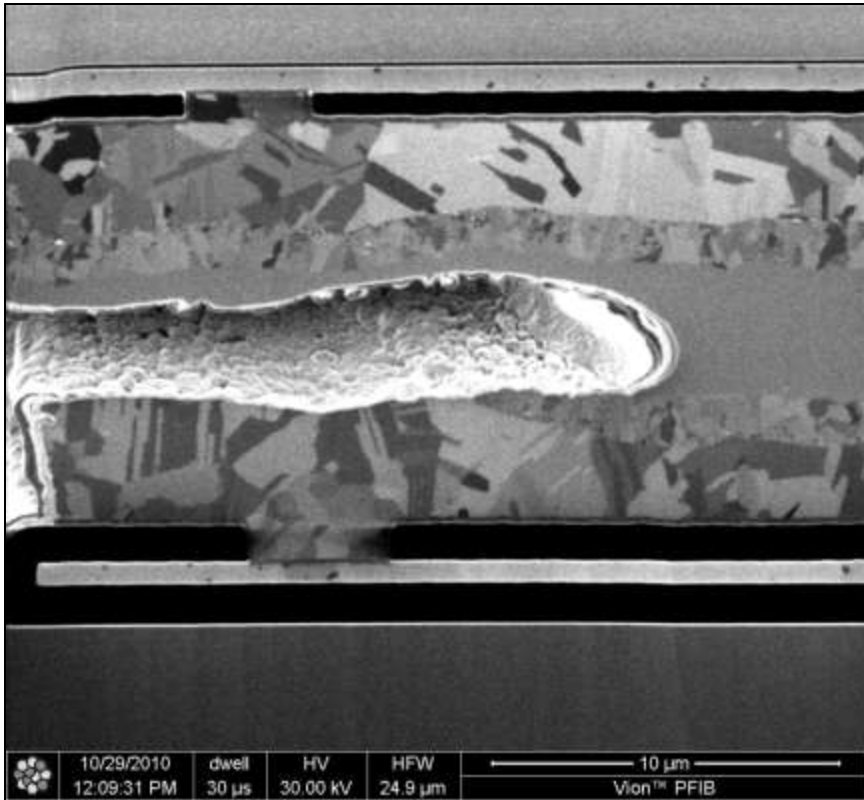


SLID bonding solid liquid interdiffusion


Reducing curtaining

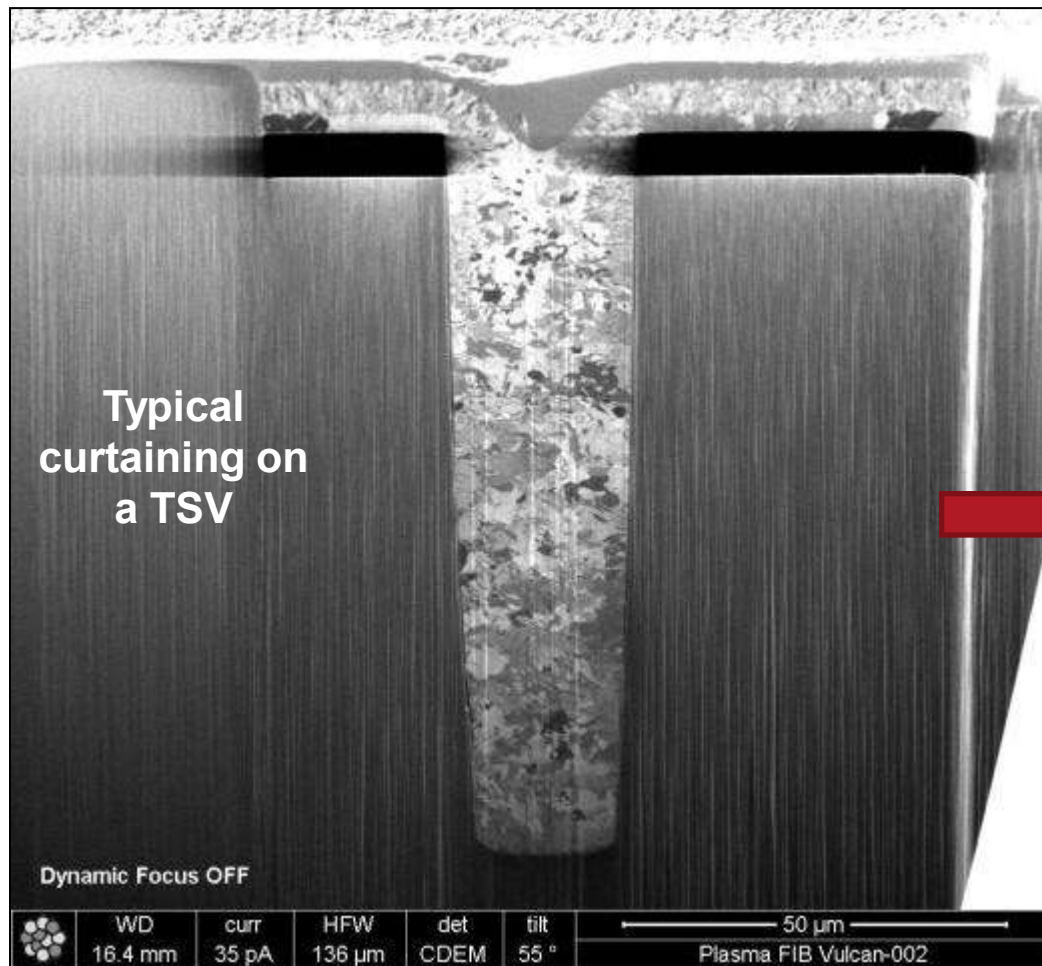

Anisotropic conductive adhesive (ACA) for Wafer-to-Wafer (W2W) bonding

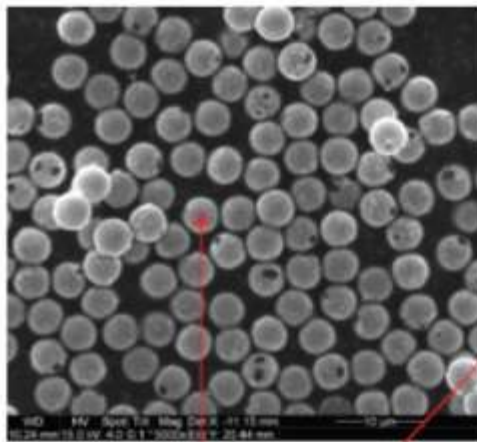


Metal coated polymer spheres (4 μm diameter)

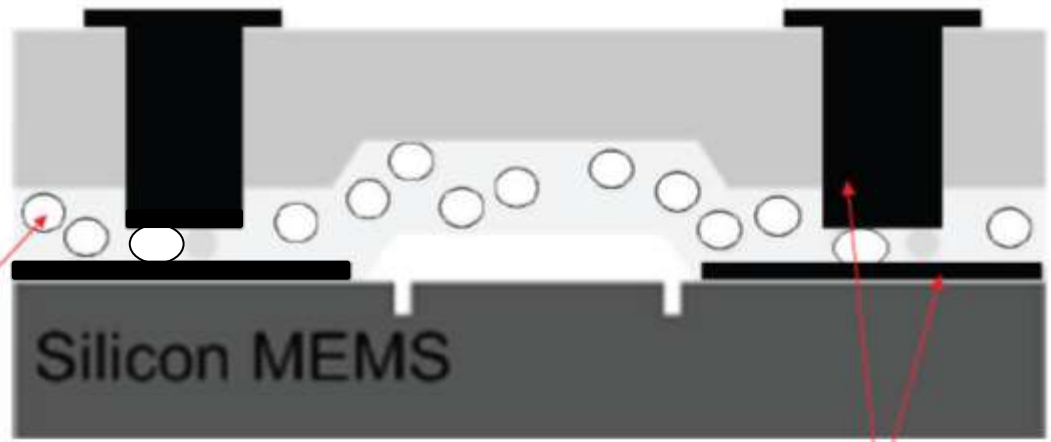



Electrical contacts

Ref: M.M.V. Taklo, T. Bakke, H.R. Tofteberg, L.G.W. Tvedt and H. Kristiansen, Proc. IMAPS Device Packaging Conf., Scottsdale, Arizona, 2011

Anisotropic conductive adhesive (ACA) for W2W bonding

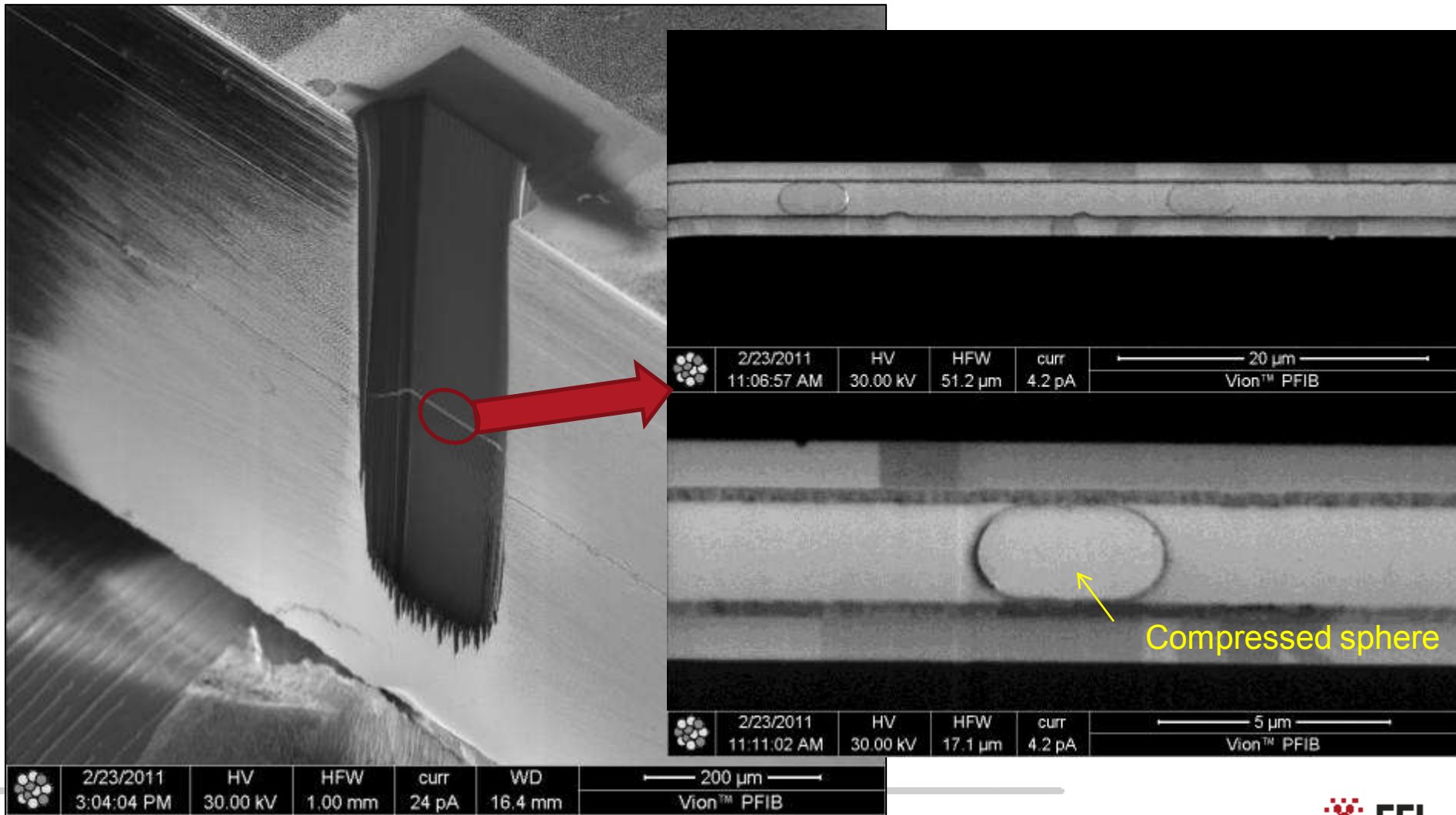

Fiber reinforced polymer bar used for tensile testing

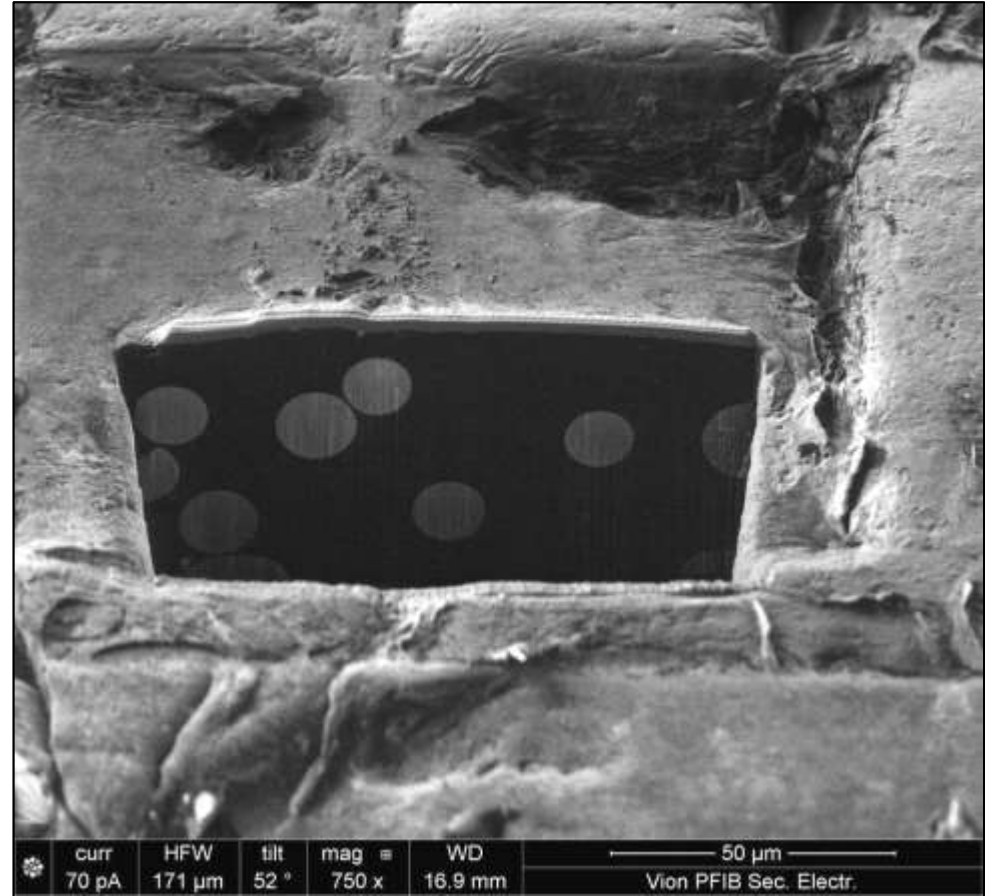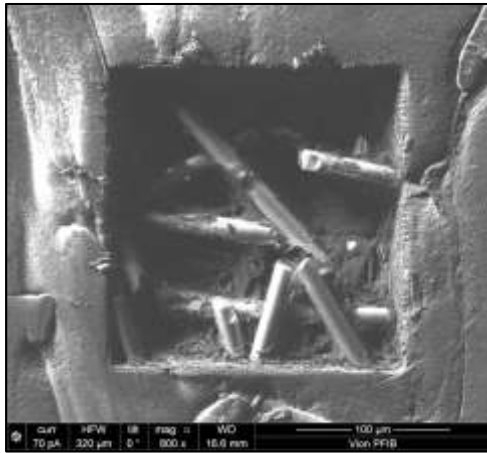

Summary

3D IC technology needs metrology and root cause analysis down to the sub-micron level for development and failure analysis

- But it takes too long with traditional Ga-FIB

Plasma FIB technology brings site specific advantages of Ga FIB to chip/package scaled problems

- More than 20x faster than traditional FIB
 - Capable of high-precision final cuts and high-resolution (sub-30 nm) imaging
 - Provides faster development feedback and failure analysis
- Non-semiconductor applications
- Many possible applications - just starting to scratch the surface
 - There are many interesting beam-sample and beam-chemistry interactions to discover
 - Ga has a 20-year head start, but PFIB learning curve on new samples is rapid!!



Acknowledgements

- FEI Engineering Team for development and commercialization of the Vion PFIB source, column and system
- Applications results and discussions: Sean Kellogg, Tom Miller, Ruud Schampers, Theresa Graupera, Chad Rue, German Franz, Kimball Skinner, Trevan Landin, Rob Routh, Remco Geurts, Loek Kwakman, Pete Carleson, Kenny Mani
- A part of the work has been performed in the project JEMSiP_3D, which is funded by the Public Authorities in France, Germany, Hungary, The Netherlands, Norway and Sweden, as well as by the ENIAC Joint Undertaking

References (1/2)

1. N.S. Smith, W.P. Skoczylas, S.M. Kellogg, D.E. Kinion, P.P. Tesch, O. Sutherland, A. Aanesland, and R.W. Boswell: “High Brightness Inductively Coupled Plasma Source for High Current Focused Ion Beam Applications”, J. Vac. Sci. Technol., B24(6), 2902-2906, 2006
2. S. M. Kellogg, A. Graupera, R. Hoelle, T. Miller, S. Zhang, D. Laur, and A. Dirriwachter: “A System for Massive, Rapid Material Removal for Device Analysis in Monolithic 3D Integrated Circuits”, Presented at 53rd International Conference on Electron, Ion and Photon Beam Technology and Nanofabrication, Florida, May 2009
3. S.M. Kellogg, R. Schampers, S.Y. Zhang, A.A. Graupera, T. Miller, W.D. Laur, and A.B. Dirriwachter, “High Throughput Sample Preparation and Analysis using an Inductively Coupled Plasma (ICP) Focused Ion Beam Source”, Microsc. Microanal., 16 (Suppl 2), 222-223, 2010
4. R.J. Young, C. Rue, M. Schmidt, R. Schampers and D. Wall: “Site-specific analysis of advanced packaging enabled by focused ion beams (FIBs)”, International Wafer Level Packaging Conference Proceedings, Santa Clara, California, October 11-14, 2010
5. R.J. Young: “Site-Specific Analysis of Advanced Packaging Enabled by Focused Ion Beams”, Electronic Device Failure Analysis, Vol. 13(1), 12-19, 2011

References (2/2)

6. P. Ramm, A. Klumpp, G. Franz, and L. Kwakman: “Failure Analysis and Reliability of 3D Integrated Systems”, Proc. IMAPS Device Packaging Conf., Scottsdale, Arizona, 2011
7. M.M.V. Taklo, T. Bakke, H.R. Tofteberg, L.G.W. Tvedt and H. Kristiansen: “Anisotropic Conductive Adhesive for W2W Bonding”, Proc. IMAPS Device Packaging Conf., Scottsdale, Arizona, 2011
8. R.J. Young: “Site-specific Analysis for Failure Analysis and Reliability of 3D Systems”, Chip Scale Review, May/June, p27, 2011
9. L. Kwakman, G. Franz, M.M.V. Taklo, A. Klumpp and P. Ramm: “Characterization and Failure Analysis of 3D Integrated Systems using a Novel Plasma-FIB System”, NIST: Frontiers of Characterization and Metrology for Nanoelectronics, Grenoble, May 2011
10. R.J. Young, C. Rue, R. Routh, G. Franz, L.F.Tz. Kwakman, P. Ramm, A. Klumpp, and M.M.V Taklo: “Site-specific Failure Analysis and Reliability Testing of 3D Systems using “Plasma” FIB”, SEMATECH 3D Interconnect Metrology Workshop, San Francisco, July 13th 2011
11. R. Young, C. Rue, S. Randolph, C. Chandler, G. Franz, R. Schampers, A. Klumpp, and L. Kwakman: “A Comparison of Xenon Plasma FIB Technology with Conventional Gallium LMIS FIB: Imaging, Milling, and Gas-Assisted Applications”, presented at M&M 2011
12. <http://www.fei.com/vion>